

# A Systems Biology Approach; Proteomic Study on Rapamycin Treated Human Cancer Cell Line

ASMS 2010

Sudha Rajagopalan, Nilanjan Guha,  
Yugandhar Reddy, Siji Joseph and Syed  
Lateef, Agilent Technologies India Pvt.  
Ltd, Bangalore India

# A Systems Biology Approach; Proteomic Study on Rapamycin Treated Human Cancer Cell Line


## Introduction

Rapamycin is an immunosuppressant drug which is recently shown to inhibit cancer growth. Rapamycin specifically inhibits mTOR activity and cellular hyper-proliferation in many cells types resulting in G1 growth arrest. A comprehensive study at genome, proteome and metabolome level will enable better understanding of the biological processes affected by this drug. Multiple instrument platforms and advanced software are required for such an analysis. In the current study, we have used systems approach to understand the effect of rapamycin treatment on various biological pathways using state of the art platforms from Agilent. Effect of rapamycin treatment was studied on human cancer cell line HEK293.

## Study Design

Human cell line HEK293 was treated with rapamycin and cells were harvested at 15min, 45min, 2hrs and 16hrs after rapamycin treatment. Control samples treated with vehicle alone were also harvested at each time point. RNA, Protein and metabolites extracted from treated and control samples were compared. Five biological replicates were performed on five different days.

Figure1 Study Design

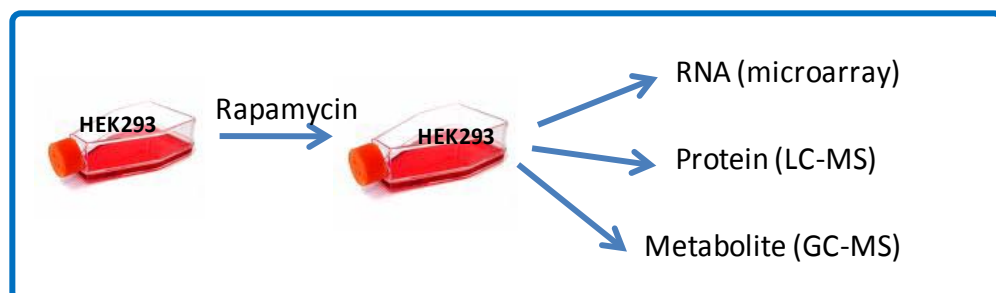

## Experimental

### Cell Culture and harvesting cells

HEK293 cells (ATCC: CRL-1573) were cultured in DMEM medium containing 10% fetal bovine serum, 10 U of penicillin/ml, and 10 µg of streptomycin/ml. Cells were grown to an approximate density of 10<sup>5</sup>/ml, 40 nM rapamycin or DMSO vehicle was added, Cells were harvested at 15, 45, 120 min and 16 h and aliquoted so that each tube contains 10<sup>6</sup> cells.

### RNA Analysis

RNA was extracted from the harvested cells using Stratagene Absolutely RNA kit (Cat# 400814). Both miRNA and mRNA species are quantitatively recovered from the Stratagene kit, allowing gene expression studies of both miRNA and mRNA from the same sample.

The RNA integrity number (RIN) values were determined for the extracted RNA using RNA 6000 Nano Lab Chip kit on the Agilent 2100 Bioanalyzer and the quality of RNA was found to be good.


25 ng of RNA samples for each control and treated samples were labeled with Cy3 using Agilent's Low Input Quick Amp kit and gene expression microarray analysis was performed using Agilent Human GE 4X44K V2 Microarray kit (G4845A). The data was analyzed using GeneSpring 11.0.2. (Agilent) Entities having a p-value ≤ 0.1 and fold change ≥ 1.4 were considered differentially expressed.

## Experimental

### Protein Analysis

Protein was extracted using PPS silent surfactant (Stratagene), reduced, alkylated and digested. Samples were analysed in five replicates on HPLC chip-QTOF system. Differential features obtained in the comparison of control and treated samples were then identified using targeted MS/MS analysis. Workflow used for the proteomics analysis is shown below.

Figure2 Workflow used for proteomics Analysis

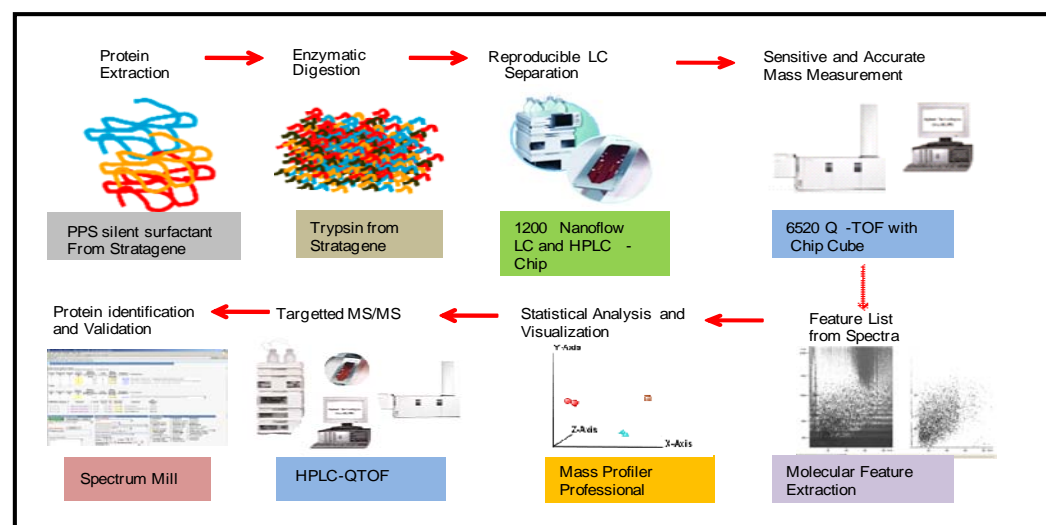

All LC/MS experiments were performed on an Agilent 1200 Series HPLC-Chip/MS interfaced to an Agilent 6520 Accurate-Mass Q-TOF LC/MS. The LC system consisted of a capillary pump for sample loading, a nanoflow pump and a thermostated microwell-plate autosampler. The HPLC-Chip configuration consisted of a 160 nL enrichment column and a 150 mm x 75 µm analytical column (Zorbax 300SB-C18). Mobile phases employed were: A). 0.1% formic acid in water and B). 90% acetonitrile with 0.1% formic acid. LC and MS conditions used were as reported earlier.<sup>1</sup>

Mass profiler professional (MPP), Agilent version B02.00 was used for statistical comparison of the LC-MS data and SpectrumMill version B03.00 (Agilent) was used for data base searching of the MS/MS data.

### Metabolome Analysis

Metabolites were extracted from one million cells with 2:1:2 methanol-water/chloroform mixture and derivatized in two steps before GC-MS analysis. First, using methoxyamine hydrochloride in pyridine to derivatize carbonyl groups, then using N-methyl- N trimethylsilyltrifluoroacetamide containing 1% trimethylchlorosilane to derivatize carboxyl, hydroxyl, amino, imino or sulphonyl groups. d27-myristic acid, was used as internal standard for the GC-MS analysis.

Agilent 7890 GC coupled with Agilent 5975C inert MSD was used for the data collection. A 30m long Agilent DB-5MS column with 10m Duragard precolumn (Part Number: 122-5532G) with 0.25 µm film thickness and 250µm diameter was used for the separation.. All GC-MS experiments were performed as described in the user manual (G1676-80000)<sup>2</sup>.

Statistical comparison between control and treated samples were performed using MPP (Agilent) Differentially expressed metabolites were identified using Agilent Fiehn GC/MS Metabolomics RTL Library (version June 2008).

## Results and Discussion

### Protein Analysis

Approximately Forty proteins were significantly (*P* value less than 0.05 and fold change greater than 10) differentially expressed between the treated and control samples at 16h of rapamycin treatment.

Table 1 A partial list of Differentially expressed peptides and the proteins from which they are derived

Peptide Sequence	Regulation	RT	m/z	Mass	Prot.Accession No.	Protein Name
(K)VDFPQDQLTALTGR(I)	up	31.79	780.9076	1560.802	6648067	Malate dehydrogenase, mitochondrial precursor
(K)AFYPEEISSMVLTK(L)	up	37.1	807.9243	1614.808	4529894	HSP70-HOM
(K)AFYPEEISSMVLTK(L)	up	37.1	538.9521	1614.808	4529894	HSP70-HOM
(K)LQMEAPHIIVGTPGR(V)	down	23.83	540.2879	1618.873	4503529	eukaryotic translation initiation factor 4A isoform 1
(K)DPVQEAWAEDVDLR(V)	down	30.14	821.8985	1642.771	33286418	pyruvate kinase
(K)DLADELALVDVIEDK(L)	up	48.23	829.4246	1657.853	5031857	lactate dehydrogenase A
(K)NQVALNPQNTVFDAK(R)	up	19.45	829.9158	1658.85	167466173	heat shock 70kDa protein 1B
(K)ILLANFLAQTEALMR(G)	down	51.94	852.4857	1703.951	18201905	glucose phosphate isomerase
(K)ILLANFLAQTEALMR(G)	down	51.94	568.6612	1703.951	18201905	glucose phosphate isomerase
(K)VAVLGASGGIGQPLSLLK(N)	up	48.08	897.0407	1793.09	6648067	Malate dehydrogenase, mitochondrial precursor
(R)AMADPEVQQIMSDPAMR(L)	up	28.4	945.4072	1889.855	5803181	stress-induced-phosphoprotein 1 (Hsp70/Hsp90-organizing protein)
(R)IFANTPDSGCVLGMKR(A)	down	32.95	641.329	1921.974	4505749	phosphofructokinase, muscle
(R)LTLYDIAHTPGVAADLSHIETK(A)	down	38.47	592.0709	2365.24	6648067	Malate dehydrogenase, mitochondrial precursor
(R)AEGSDVANAVLDGADCIMLSGETAK(G)	down	42.8	1247.544	2494.144	33286418	pyruvate kinase
(R)AEGSDVANAVLDGADCIMLSGETAK(G)	down	42.83	832.0589	2494.144	33286418	pyruvate kinase
(R)VETGVLKPGMVVTFAPVNVTTTEVK(S)	down	38	629.5953	2515.384	4503471	eukaryotic translation elongation factor 1 alpha 1
(K)SGDAIVDMVPGKPMCVESFSDYPPPLGR(F)	down	41.35	999.131	2995.4	4503471	eukaryotic translation elongation factor 1 alpha 1
(K)EIAEAYLGYPTNAVITVPAYFNDSQR(Q)	up	48.17	1001.162	3001.494	167466173	heat shock 70kDa protein 1B

Figure 3 PCA plot on replicate LC-MS control and rapamycin treated (16h) samples, reveals that the variance in the data can be explained by which class it came from

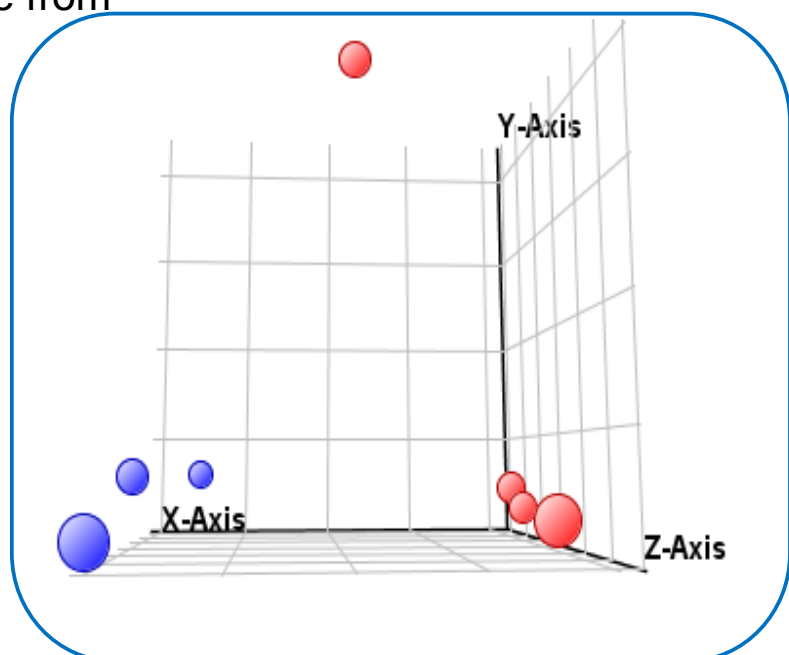


Figure 4 Mass spectrum between 29.2-29.4 minutes showing one of the differentially expressed peptides DPVQEAWAEDVDLR derived from the protein pyruvate kinase


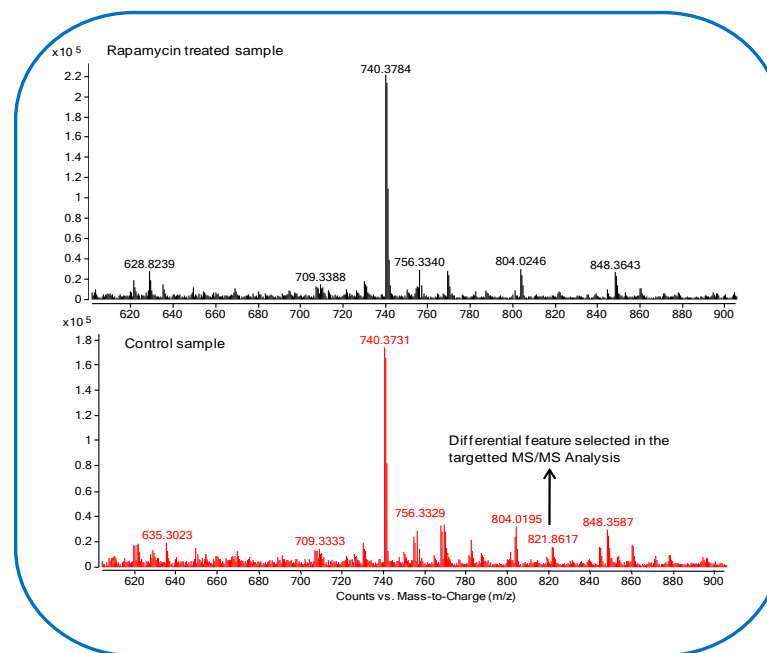


Figure 5 MS/MS spectrum of the peptide DPVQEAWAEDVDLR acquired by targeted MS/MS

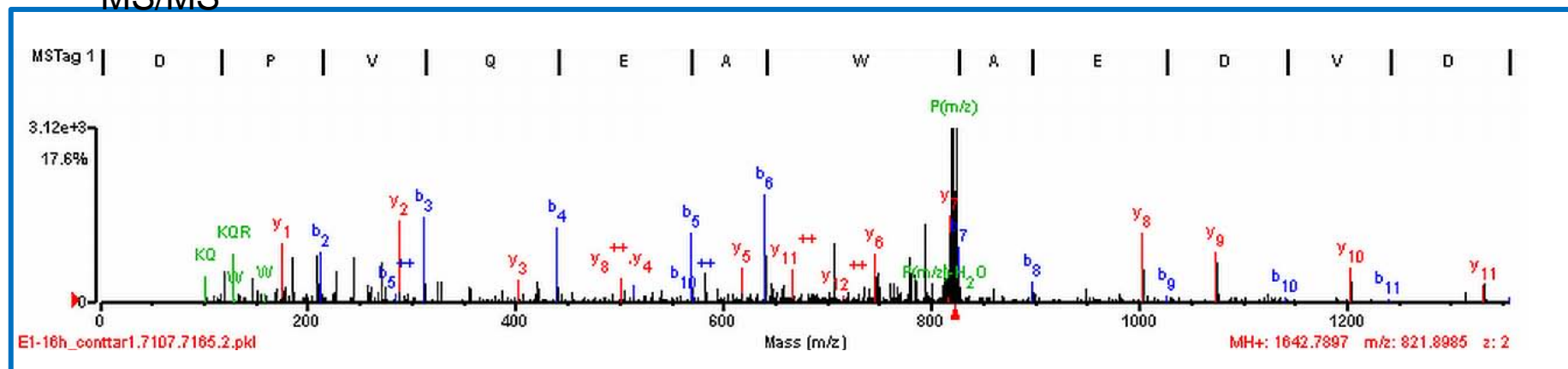

The peptide DPVQEAWAEDVDLR is co-eluting with several other high abundant peptides. This peptide may not get selected in the data dependant MS/MS analysis.

## Results and Discussion

### Genome Analysis

Approximately 72 genes were differentially expressed between the control and rapamycin treated samples at 16h. Gene expression analysis suggests that transcription, translation, cell growth, and the stress response are some of the biological processes affected by rapamycin treatment.

Table 2 A Partial list of differentially expressed genes

Fold change	Regulation	Gene Symbol	Protein name
1.6216595	up	ATF3	activating transcription factor 3
1.4631783	down	C11orf52	chromosome 11 open reading frame 52
1.6093587	down	TMEM86A	transmembrane protein 86A
1.5000856	up	HSPA9	heat shock 70kDa protein
1.434173	up	ZFP36L1	zinc finger protein 36, C3H type-like 1
1.5229359	up	HSPH1	heat shock 105kDa/110kDa protein 1
1.5675275	up	INPP5F	inositol polyphosphate-5-phosphatase F
1.5505759	up	ETNK1	ethanolamine kinase 1
1.4378525	up	TBCE	tubulin folding cofactor E
1.5906347	up	GNPDA2	glucosamine-6-phosphate deaminase 2
1.4255726	up	MRPL19	mitochondrial ribosomal protein L19

### Metabolome Analysis

Approximately 30 metabolites were found to be differentially expressed between the control and rapamycin samples at 16h. The metabolite list included up-regulation of several sugars and amino acids.

Table 3 A Partial list of differentially expressed metabolites

Compound name	[r] Vs [c] regulation	Fold change	Mass	Retention Time
gluconic acid lactone 1	down	16	147	17.66
malonic acid 1	down	16	73	8.45
D-glucose 1	up	16	186	17.63
talose 1	up	16	205	17.66
N-acetyl-D-glucosamine 1	up	16	147	18.743
10-hydroxydecanoic acid	up	16	227	15.71
glycerol 1-phosphate	up	16	299	16.04
Creatinine	up	16	115	13.68
fructose 1	up	16	133	17.17
D (+)altrose 1	up	4.92	189	17.7
L-threonine 2	up	3.98	73	11.54
L-pyroglutamic acid	up	3.96	156	13.24
glycine	up	3.79	174	10.46
2-amino-1-phenylethanol	up	3.27	174	9.42
L-(+) lactic acid	up	3.27	73	6.8

Figure 6: Total ion chromatogram from 17.6 - 19.5 minutes showing differential expression of one of the glucose derivatives

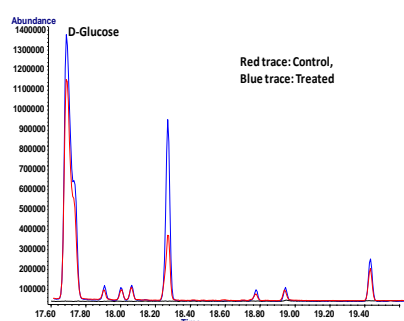


Figure 7 Fragmentation pattern observed in the experiment matches well with the library, confirming glucose identification.

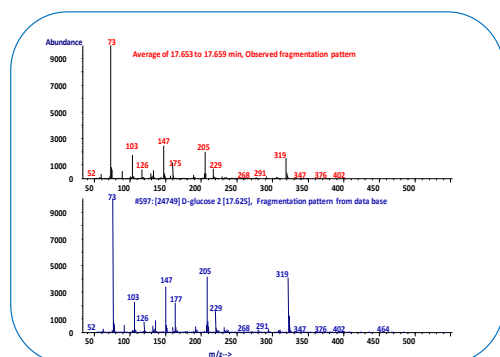


Table 4 Biological processes affected by rapamycin

Biological Process	RNA (microarray)	Protein (LC-MS)	Metabolite (GC-MS)
Glycolysis/ gluconeogenes	• Glucosamin e-6-phosphate deaminase 2	• Pyruvate kinase • Glucose phosphate isomerase • Malate dehydrogenase,	• Glucose • Fructose
Stress Response	Heat shock 70kDa protein	Heat shock 70kDa protein	----
Protein Synthesis	----	Eukaryotic translation initiation factor 4A	Glycine, Cysteine, Threonine, Valine, Leucine, Alanine

Green : Elevated                      Red: Suppressed

## Conclusions

- Glycolysis/gluconeogenesis pathway is reported to be affected by rapamycin treatment<sup>3</sup>. In this study, down regulation of pyruvate kinase, glucose phosphate isomerase and up-regulation of malate dehydrogenase are observed, indicating glycolysis is reduced and gluconeogenesis is enhanced due to rapamycin.
- Up-regulation of the gene for glucosamine-6-phosphate deaminase 2 and elevation of glucose in metabolite analysis supports this observation.
- Up-regulation of Heat shock proteins, which is known to be a stress response is observed in both proteomics and gene expression analysis.
- Rapamycin treatment is reported to result in decreased protein levels<sup>4</sup>. Observation of down regulation of eukaryotic translation initiation factor 4A in proteomic analysis and up-regulation of amino acids in the metabolite analysis indicates down regulation in protein synthesis.
- Proteomic, genomic and metabolomic analysis using various platforms enables better understanding of biological process.

### Future Direction

Preliminary results obtained in this study is encouraging and further experiments on more biological replicates are currently in progress. Integration of data from genome, proteome and metabolome analysis to elucidate biological processes will be available in future versions of Agilent GeneSpring software.

### References

- Agilent Publication 5990-4404EN Profiling Approach for Biomarker Discovery using an Agilent HPLC-Chip Coupled with an Accurate-Mass Q-TOF LC/MS.
- Agilent Publication G1676-90000 Agilent G1676AA Fiehn GC/MS metabolomics RTL library user guide.
- Houde VP, Brûlé S, Festuccia WT, Blanchard PG, Bellmann K, Deshaies Y, Marette A. Diabetes. 2010 Mar 18.
- Tao Peng, Golub, TR and Sabatini DM Mol.Cell.Biol. 2002, p. 5575-5584.