

## Base Neutral 4

### Column: DB-5.625

30 m x 0.25 mm I.D., 0.25 µm

J&W P/N: 122-5631

Carrier: Helium at 33 cm/sec

Oven: 35°C for 5 min

35-280°C at 8°/min

280°C for 20 min

Injector: Splitless, 250°C

Detector: FID, 300°C

Nitrogen makeup

gas at 30 mL/min

### Compounds and Retention Times:

39. 4-Chlorophenylphenyl ether (24.85)

40. Benzidine (32.11)


41. Di-n-octyl phthalate (38.77)

42. Benzo(k)fluoranthene (39.98)

43. Benzo(a)pyrene (41.26)

44. Indeno(1,2,3-c,d)pyrene (47.82)

45. Benzo(g,h,i)perylene (49.67)

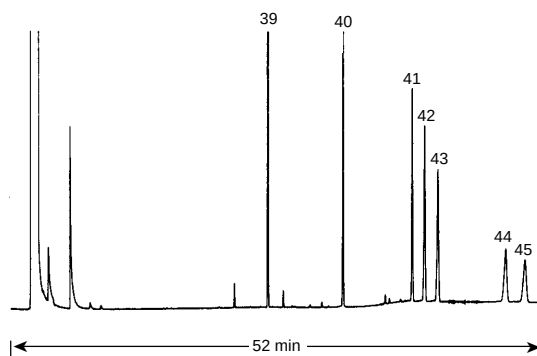

C021