

Sedative Hypnotics

Column: DB-5ms EVDX
12 m x 0.20 mm I.D., 0.33 μ m

J&W P/N: 128-8512

Carrier: Helium at 35 cm/sec, measured at 55°C

Oven: 55°C for 1 min

55-130°C at 25°/min

130-325°C at 15°/min

325°C for 4 min

Injector: Splitless, 250°C


45 sec purge activation time

1 μ L of 50-100 ng/ μ L standard in methanol

Detector: MSD, 280°C transfer line

full scan at m/z 35-400

1. Ethchlorvynol
2. Ethinamate
3. Pyrrithyldione
4. Talbutal
5. Meprobamate
6. Glutethimide
7. Methaqualone
8. Propiomazine
9. Haloperidol
10. Sulpiride
11. Droperidol

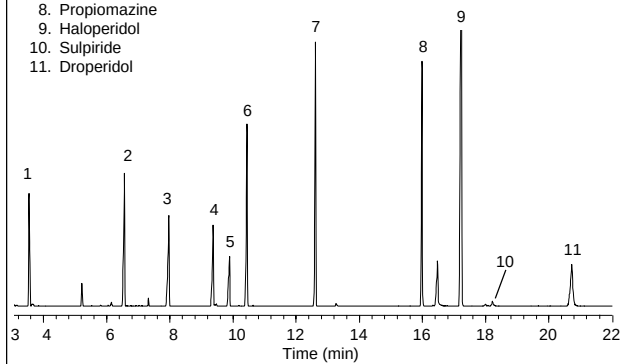

C218